HUMAN CARE
Your Life. Your Way.

Bedtable

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User Manual

Important!

You must read the user manual for your device prior to usage. Keep this booklet and information handy for future use.

Symbols

WARNING symbol indicates a potentially hazardous situation which, if not avoided, could result in serious injury to the user or damage to property and/or the device itself.

INFORMATION symbol indicates recommendations and information for proper, trouble-free usage.

WEIGHT CAPACITY symbol indicates maximum user weight for the product. Do not exceed this weight in any circumstance.

CE MARKING symbol indicates product meets requirements of EU Directive 93/42/EEC (MDD) and/or EU Regulation 2017/745 (MDR).

MEDICAL DEVICE symbol indicates that the product is a medical device according to the definition of medical device in EU Directive 93/42/EEC (MDD) and/or EU Regulation 2017/745 (MDR).

RECYCLE according to local regulations.

READ INSTRUCTIONS and ensure you fully understand them before using this product.

Product modification is not permitted.

Before use, always check that the sling is not damaged or worn.

Human Care's products are continuously being developed and updated to provide our customers with the highest quality. We reserve the right to make design changes without prior notice.

Always make sure that you have the most recent version of the manual which is available for downloading from our website at www.humancaregroup.com.

Contact your local distributor if you have any questions about the product and its use.

Notice to user/patient in case of serious incident

Any serious incident that has occurred in relation to the product, should be reported to the local contact, who reports to the manufacturer, and the competent authority of the country in which the user/patient is established.

Compliances and Standards

Human Care is an ISO 13485:2016 certified Swedish medical device company. The Quality Management System is in compliance with US 21 CFR part 820.

The product is CE marked in accordance with EU

Regulation 2017/745 (MDR), as class I medical device.

The CE mark is on the product.

About our products

On delivery

Check the following before using your new product:

1. Check that no damage has occurred during transport (if this is the case, please immediately inform the delivery agent).
2. Check that the number and type (s) of products detailed on the delivery certificate agree with the goods delivered.
3. Check that the delivery is of the desired model and type of product.

Please contact us immediately in case of any delivery problem.

Return or refund

If any product has to be returned for any reason, please contact us to obtain a claim number and advice on how to proceed.

The goods should be returned in complete and undamaged condition in their original packaging.

The following information should accompany any product on its return. The product's identification number, delivery date, the date the user took possession of the aid, the date of the product's last service, the person whom we should contact about the fault, whether it be the user or another, and the nature of the fault.

Human Care HC AB reserves the right to inspect the returned product. If no fault can be found, the item will be returned to the user. Transport costs may then be chargeable.

Product lifetime

The product's expected lifetime is seven (7) years, if the product is used as intended and maintained according to the manufacturer's instructions, depending on the intensity of use and maximum load applied during use. If the product label is no longer legible, the product should be discarded.

Bedtable

Product label

The approved product label is attached to the table. Before use and after reconditioning, make sure the label is present and legible.

The label is located on the foot of the table.

	Article number
	Maximum user weight
	Date of manufacture
	Address of manufacture
	Unique Device Identification, including (01) EAN code, (11) production date, (21) serial number
	Serial number

	Read the manual
	Product width
	Indoor / outdoor use
	CE marking
	MD Medical Devices

Additional Warnings

- Use the table only for its intended purpose. Do not use it to support a person, for transfers, or as a walking aid.
- Do not remove or bypass safety components. Inspect the controls and castors at least weekly. Before first use, read this manual and make sure that the user or carer understands all controls.
- Use the table only on a firm, level floor. Keep the movement area clear of obstacles.
- Lock the castors before the user leans on, eats from or reads at the table.
- Do not place the product next to open flames, heaters or other heat sources.
- Keep fingers, clothing and bedding clear of the height and tilt mechanisms when adjusting the table.
- Maintenance and repairs must be carried out by competent personnel using approved parts.
- Do not exceed the 20 kg safe working load. Place loads centrally and distribute them evenly.

Intended Use

The K058PTS and K058PTS-C are mobile, height-adjustable and tiltable patient tables for hospitals, care facilities and home-care environments.

Use the product indoors at temperatures from 0 to +40 °C. It is intended for users aged 12 years and older, with assistance when required by the user's condition or local care plan.

Use only the accessories supplied with the table or accessories approved by Human Care. The maximum permitted load is 20 kg, evenly distributed over the tabletop

Read this user manual before using the product.

Safety devices

Before each use, check that the table is stable, the tabletop is secure and the adjustment controls operate correctly.

1. Position the table beside the bed or chair and make sure there is enough clearance for adjustment.
2. Press each brake pedal fully down to lock the castor. Lift the pedal to unlock it.
3. Support the tabletop with one hand. Loosen the height handgrip, move the top to the required height, then tighten the handgrip firmly.
4. Support the tabletop before operating the tilt handle. Select the required angle and make sure the mechanism is locked before placing items on the surface.
5. After adjustment, gently test the tabletop and base for movement before use.


20 kg maximum permitted load. keep the load central and evenly distributed.

General features and accessories

- Manual height adjustment with locking handgrip
- Tiltable eating and reading surface
- K058PTS: 50 mm castors, two with brakes
- K058PTS-C: 75 mm enclosed twin castors, four with brakes Hygienic ABS tabletop
- Powder-coated metal frame

Bedtable

Assembly

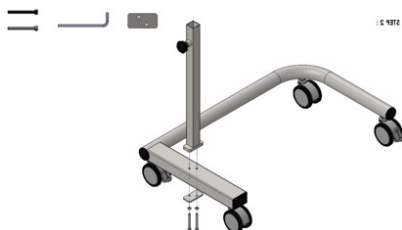
1. Place the outer column over the holes in the base. Turn the star knob outwards.

Place the mounting plate under the base.

Install the M6 spring washers and M6x55 Allen screws.

Tighten the screws alternately with the Allen key until the column is firmly seated.

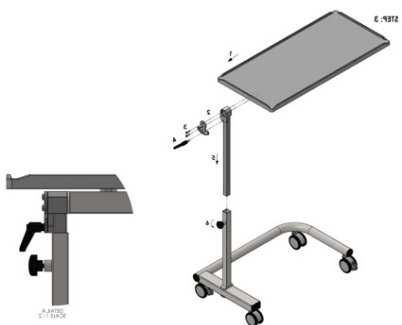2. Place the table top on the upper bracket of the inner column.

Install the coupling plate and the Allen screws. Tighten with the Allen key.

Install the tilt mechanism handle as shown in the detail drawing.

Lower the inner column into the outer column while supporting the table top.

Select the height and tighten the pre-assembled star knob.

Clamp risk: Support the tabletop during the entire assembly and keep fingers away from the column openings and tilt mechanism.

Check before use (the table must not be used if the final inspection cannot be approved):

- That the pre-assembled wheels are securely in place and that the brakes work. The outer column is flat against the base and all screws are tightened. The tabletop is secure and the tilt handle locks in the selected position.
- That the height knob locks the tabletop without play.
- That the table stands stably on a flat floor and no parts are damaged.
- That the product labeling is in place and legible.

Technical details

Specification	Value
Height range	680-1000 mm
Tabletop tilt angle	-40° to +40°
Tabletop length	820 mm
Tabletop width	380 mm
K058PTS castor diameter	50 mm
K058PTS-C castor diameter	75 mm
Maximum load	20 kg
Product weight	11.6 kg


Bedtable

Spare parts

Part	Article no.	Image ref.	Model
Eating/reading board	K058-PT/TT	1	Both
Height adjustment handgrip	K058PTS/HIGHTADJ	2	Both
Lock handle	PKC630830	7	Both
75 mm enclosed brake castor	K058-CWL	3	K058PTS-C
75 mm enclosed castor	K058-CW	4	K058PTS-C
50 mm castor without brake	K058PTS/Castor-nolock	5	K058PTS produced before 2026-10
50 mm castor with brake	K058PTS/Castorlock	6	K058PTS produced before 2026-10
50 mm castor without brake new	94432	5	K058PTS produced after 2026-10
50 mm castor with brake new	94433	6	K058PTS produced after 2026-10

Safety devices

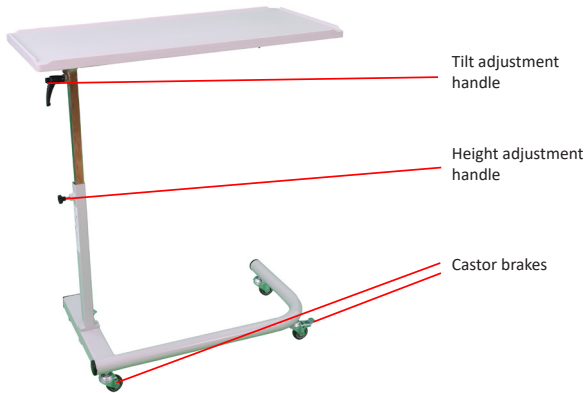Care instructions

Check for wear and tear

Check the product for wear and tear.

Do not use damaged products.

Transfer products are subject to a natural ageing process. Inspect the equipment regularly to ensure proper functionality.

Cleaning

Wipe the frame and tabletop with a damp cloth. Use commonly available cleaning and surface disinfecting agents that are compatible with ABS plastic and powder-coated metal.

Do not use abrasive pads, solvents or cleaners that can damage the surface. Rinse or wipe away cleaning residues and dry the product before use.

Before high-pressure cleaning. Lock the castors and secure the tabletop. High-pressure cleaning may only be used when permitted by the approved local reconditioning procedure.

Store/transport

Store product in temperatures between -30 to +40 degrees celsius.

Warranty and support

If you need information or support, please contact your Human Care office, www.humancare-group.com, or distributor.

Recycling

Once the product has reached the end of its life, it should be disposed according to your country's rules. The table is designed so its principal material groups can be separated at end of life. Follow local requirements for recycling medical equipment. Sort the tabletop and plastic handles as plastic where accepted. Sort the frame and metal fittings as metal. Separate castors into metal, plastic and rubber streams where local facilities allow. Sort cardboard and plastic packaging separately.

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